

BİOREDRLIGHT 360Z PATENTED TECHNOLOGY



ÜRÜN DETAYLARI

Güzellik salonlarına veya sağlık tedavi kulüplerine herhangi bir randevu almanıza gerek yok, güzellik ve zindeliğe giden yol, kendi evinizin sınırları içinde, istediğiniz zaman geçilebilir. Mavi ve yeşil de dahil olmak üzere 7 dalga boyu bulunmaktadır. Kişisel kullanım için ideal bir cihazdır. Ayrıca güzellik salonlarında ve fizik tedavi merkezlerinde bölgesel tedavi için kullanılabilir

Dalga boyu: Kanal 1 # Mavi (450nm)

Kanal 2# Yeşil (528nm)

Kanal 3 # KIRMIZI (630nm + 660nm)

Kanal 4 # NIR (850nm + 830nm + 810nm)

BİOREDRLIGHT ın gelişmiş RP serisi, yarım vücut terapisi sunarak sağlık ve zindeliğe bütünsel bir yaklaşım sağlar. 6 varsayılan mod ve 1 özelleştirilebilir mod ile sürekli ayarlama zahmetine girmeden sorunsuz bir deneyimin keyfini çıkarabilirsiniz.

Yenilikçi çift modlu dokunmatik tuşlarımız ve uzaktan kumandamız, cihazı uzaktan zahmetsizce yönetmenize olanak tanıyarak parmaklarınızın ucunda basitlik ve kontrol sağlar. Bu dönüştürücü deneyimi benimseyin ve getirdiği olağanüstü avantajların tadını çıkarın!



RP360Z

LEDs:	72x5W 3535 LED
Wavelength:	Channel 1# RED (630nm+660nm)
	Channel 2# NIR (850nm+830nm+810nm)
	Channel 3# Blue (450nm)
	Channel 4# Green (528nm)
Beam angle:	30 degrees
Control way:	Touch screen & Remote controller
Voltage:	AC85-240V
Irradiance:	>520mw/cm ²
Pulse:	0-5000Hz
Flicker:	Flicker free
EMF:	0.0μt
Use way:	Table stand
Lifespan:	>100,000 Hours
Warranty:	3 Years
Material:	Iron
Product size:	300x300x60mm / 12"x12"x2.4"
Product weight:	5KG/11LB
Gross weight:	5.5KG/12.1LB

For RP Series

Menu Introduction



Skin care

Red&NIR&Green, 20mins
(10 mins R/NIR/G + 10 mins R/NIR)



Fat burning

Red&NIR&Pulse, 30mins



Recovery

Red&NIR&Pulse, 30mins



Sleep

Red&NIR, 20mins



Inflammation

Red&NIR&Blue&Pulse, 30mins
(15 mins R/NIR/B + 15 mins R/NIR)



Health care

Red&NIR, 30mins



Customized mode

Set customized values of Red/NIR/Blue/Green/Time/Pulse and save for quick use.

Getting Started

Plug in the device and switch it on, the screen will turn on.

Default modes:

Choose your mode, then start the corresponding therapy.




Customized mode:


- Press "start" to begin the saved settings.
- Press "settings" to set your values of spectra, pulse, time, and press "done" to save.



3.5 inch high-definition touch display
6 Default modes + 3 customized modes

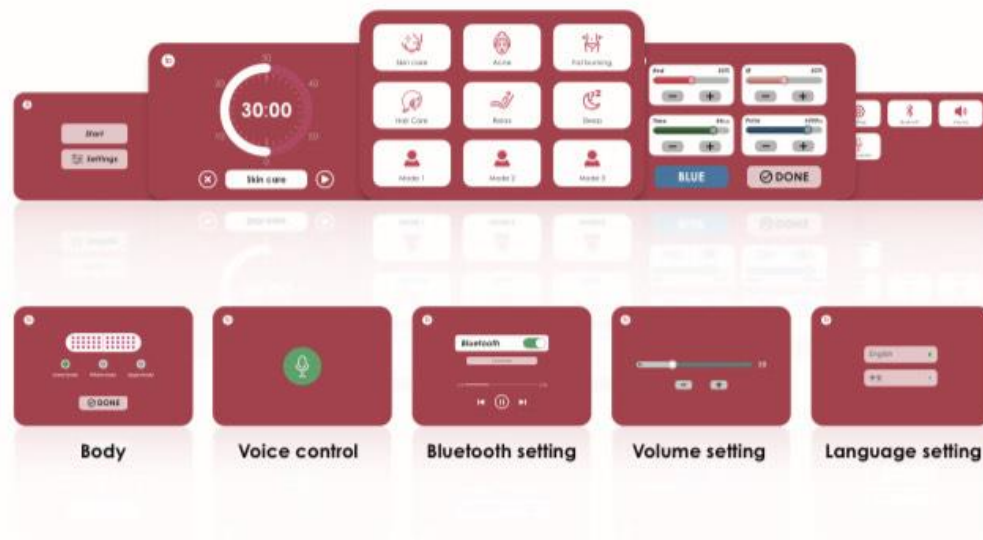

IRRADIANCE TEST

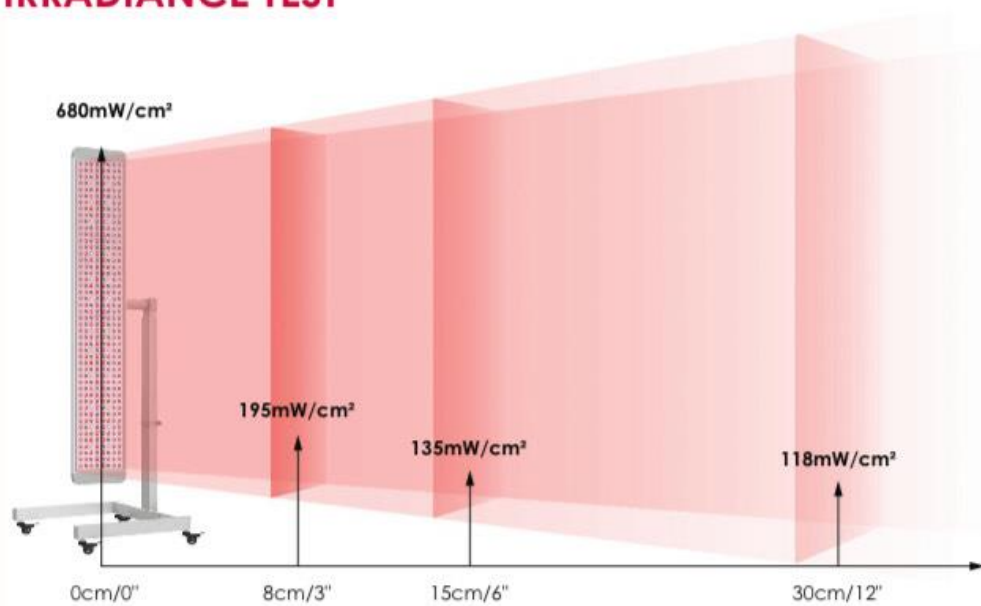


BioRedLight TOUCH SCREEN


6 Default modes + 1 Customized mode





BioRedLight REMOTE CONTROLLER

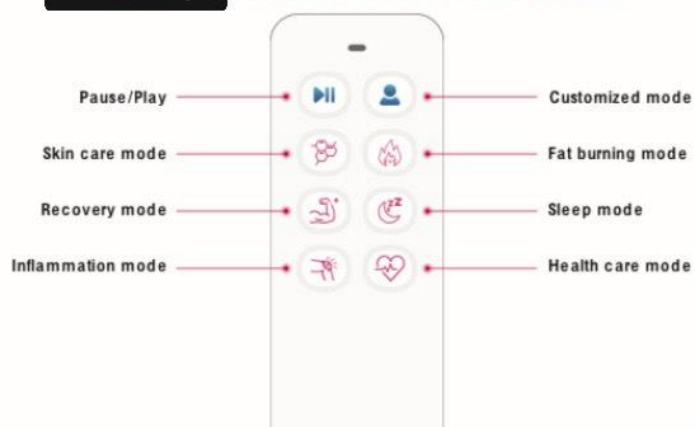




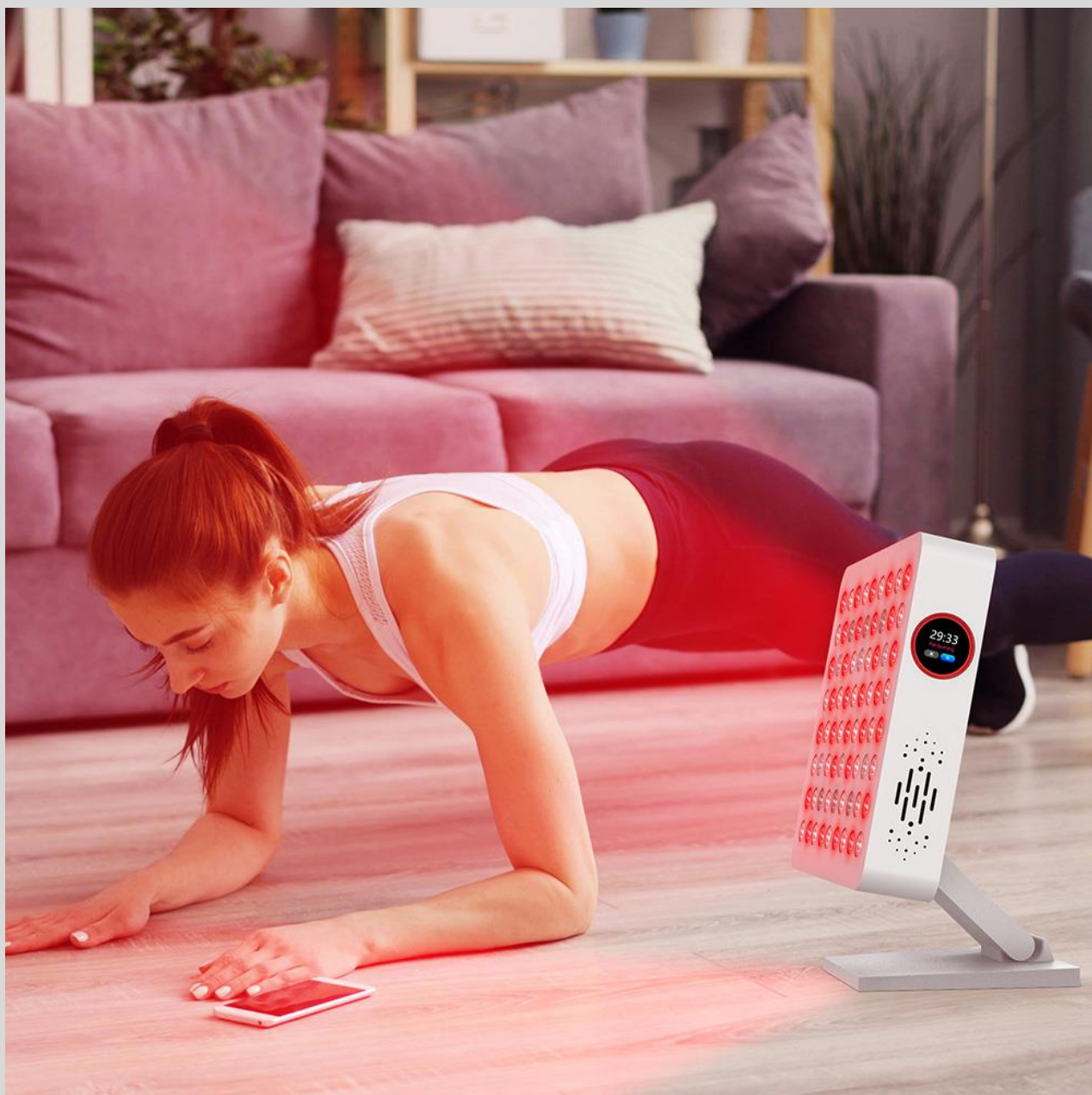

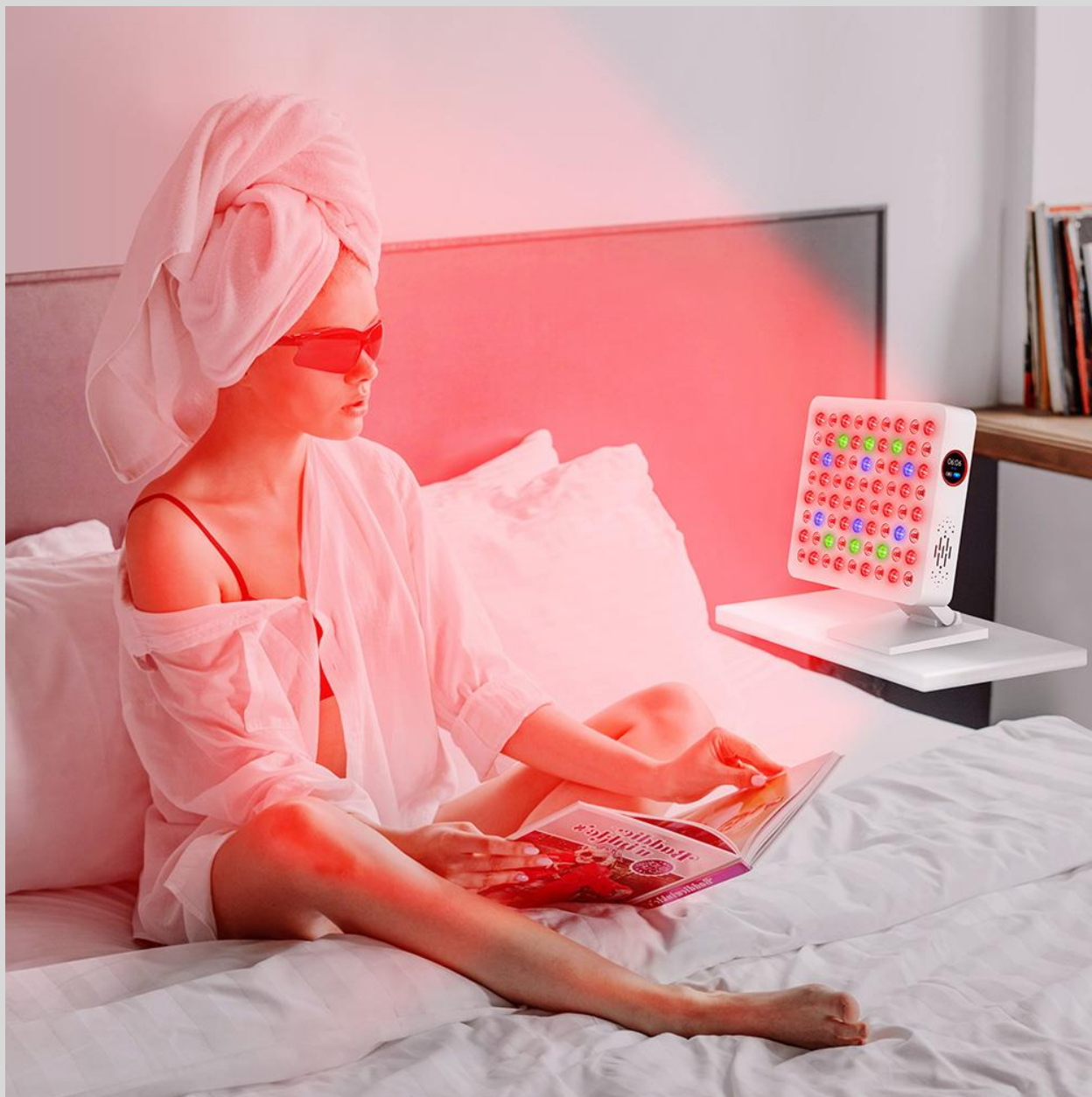

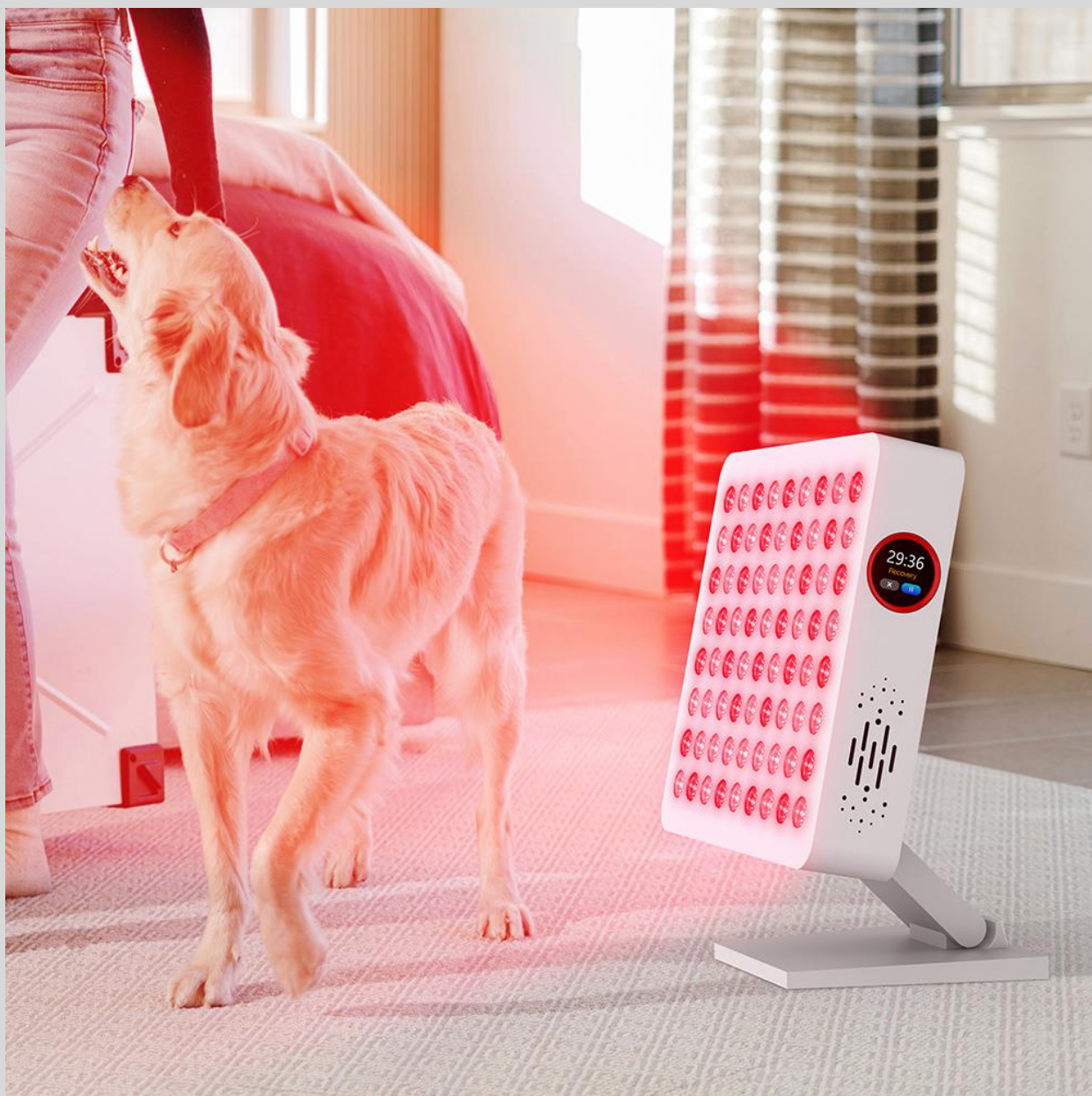

BioMol RL Series

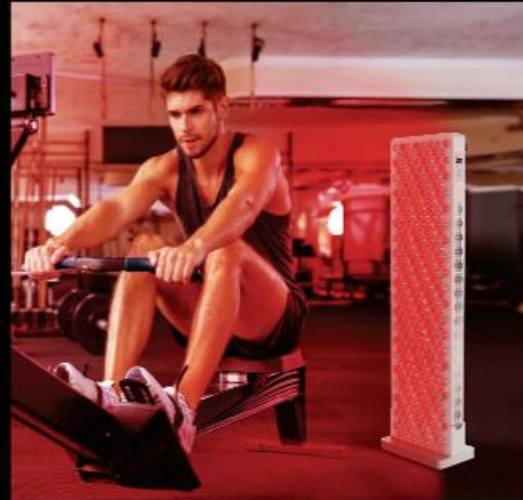

⚠ Please wear glasses or close eyes when enjoying therapy.